

Instructions for 1323 AP / AL

Polycentric 5 Bar

Pneumatic Knee Joint



ST&G USA Corp.
2691 Saturn St.
Brea, CA 92821

Phone: (714) 524-0663
Fax: (714) 364-8113
www.stngco.com

1. Description and purpose


These instructions are for use by the practitioner.

- The 1323AP/AL knee is to be used exclusively as part of a lower limb prosthesis
- Recommended for amputees with K3.
- Weight limit for a user is up to 125 kg / 275 lbs

Contra-indications

- Residual muscular weakness, contractures or proprioceptive dysfunction including poor balance.
- Contra lateral joint instabilities or pathology
- Complicated conditions involving multiple disabilities

⚠ Ensure that the user has understood any Instructions for use, drawing particular attention to the safety information.

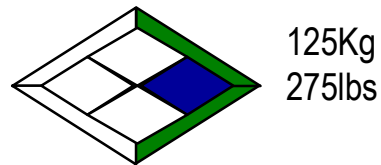

Product Code

1323AP/AL

Polycentric Pneumatic Knee Unit (Aluminum) with Stance Flexion Control

2. Construction

Principal Parts

Frame	Aluminum Alloy, Brass, Stainless Steel, Steel
Knee head	Aluminum Alloy, Stainless Steel
Knee control	Various materials principally Aluminum Alloy, Stainless Steel, Poly Urethane Copper

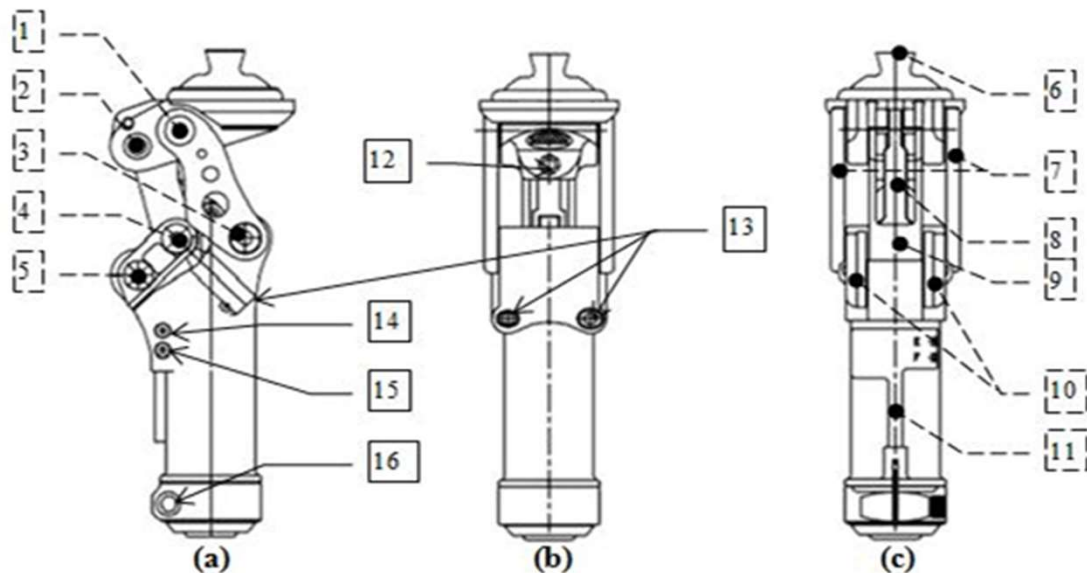


Fig. 1 (a) Side View (b) Frontal View (c) of Knee Unit

- 1) The First Axis
- 2) The Second Axis
- 3) The Third Axis
- 4) The Fifth Axis
- 5) The Fourth Axis
- 6) Knee Head
- 7) Side Bars
- 8) Piston Rod
- 9) Back Linkage
- 10) Fifth Axial Bars
- 11) Knee Body
- 12) Knee Head Level Adjusting Screw
- 13) Flexion Control Adjusting Screw
- 14) Extension Resistance Adjusting Screw
- 15) Flexion Resistance Adjusting Screw
- 16) Tube Clamp Screw

Function

- The stance flexion control angle up to 12 degree for mimicking normal knee flexion from heel strike to foot flat of a gait cycle
- Pyramid and Lotus Adapter's available
- Adjustable extension assist
- Adjustable level of knee head for various users
- The most light weight in the present market



Important:

DO NOT adjust out stance flexion to 0° - must maintain minimum of 5° stance flexion up to 12°, 0° stance flexion can potentially lead to knee failure. Excessive stance flexion can lead to unwanted noise and excessive knee stance stability which will make it difficult to initiate flexion.

4 Safety Information



The Caution symbol highlights safety information which must be followed carefully.



Be aware of finger trap hazard at all times



Any changes in performance of the knee e.g. instability or lag in transition from full stance flexion moment to full knee extension moment in the knee should be immediately reported to the Clinician / Practitioner



Any excessive changes in heel height may adversely affect the stability of the knee.



The user should be advised to contact their Clinician / Practitioner if their condition changes.

5 Maintenance

- Maintenance must be carried out by qualified personnel.
- Bi-Annual inspection to be sure knee function is satisfactory.
- Check for visual defects that may affect proper function.
- A loaner system is available should servicing be required.

The wearer should be advised:

Any changes in performance of this device must be reported to the Clinician / Practitioner.

Changes in performance may include:

- Inability to lock during weight activation
- Any unusual noises

Cleaning:

- Use a damp cloth and mild soap to clean the outside surfaces.
- DO NOT use aggressive cleaning agents.
- DO NOT use any petroleum lubrication on pivots as this will void the warranty, and could attract excessive dust or aggregate which can accelerate wear and tear on the component.
- If the limb/knee comes into contact with salt or chlorinated water, it should be rinsed with fresh water and dried.

6 Limitations on use

Intended Life:

- Service life of the product is covered by the warranty period (1 year)
- This product is recommended for use with other ST&G Products.

Lifting Loads:

Amputee weight and activity is governed by the stated limits.

Combined amputee, and carrying load, should not be at, or exceed stated weight limit.

Environment:

Store at room temperature, and avoid extremely hot, cold, or moist storage environments Storage Temperature

Range: between temperatures of -10°C to 50°C [14°F and 122°F]

Do not store in wet or very hot environment.

Allowable environmental conditions

Use – relative humidity: 0 % to 90 %, non-condensing

Storage/transportation – relative humidity: 20 % to 90 %, non-condensing

Splash - Fresh water, rain, salt water, urine, dust, sand, particles of foam cosmetics

cleaning required after contact with salt-laden air, salt water, chlorinated water, urine, sand

Unallowable environmental conditions

Storage/transportation: high vibrations, impacts

Cleaning agents containing solvents



7 Alignment and Set-Up

! Users be aware of potential finger trap hazard


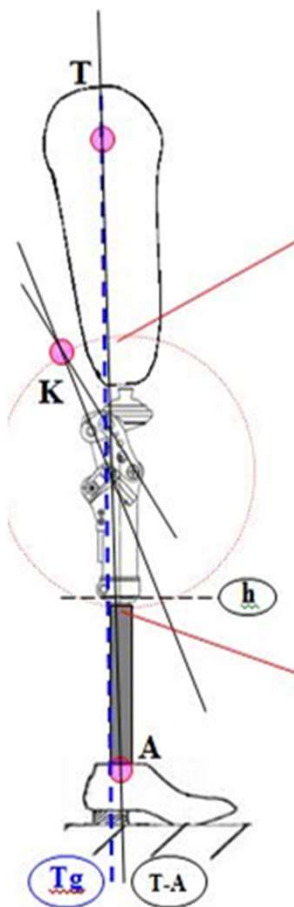


Fig. 2

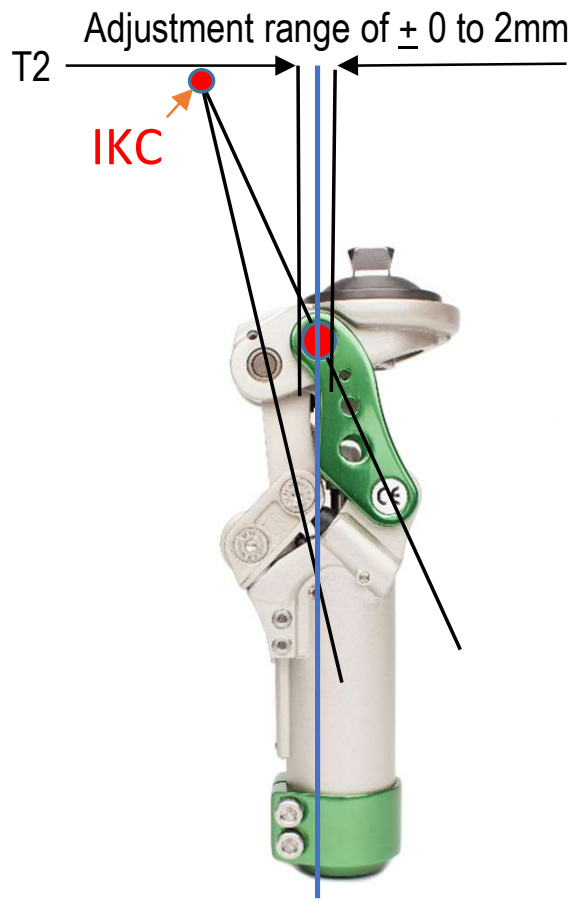

Tg
Fig. 3



Fig. 4

! Note: 4-bar knees inherently are very stable due to the geometry built into each design. This is commonly referred to as the Instant Knee Center (IKC). The IKC point when doing bench alignment, will fall behind the traditional TKA line that we will reference. (Fig. 2,3) Tg line in Fig. 3 is ideal placement, but in certain instances, it may be necessary to accommodate placement anteriorly (into the T1 zone up to 2mm), or posteriorly (into the T2 zone up to 2mm). The Tg line is referencing a moving weight bearing line, so it could be in T1, neutral, or in T2 zones.

BENCH ALIGNMENT:

a) With prosthesis assembled, taking into account hip flexion contractures, abduction, Line Of Progression, and toe out (Fig.2), the TKA plumb line should pass through the knee center (at the proximal/anterior pivot – red dot pivot Fig.3) and in front of the K point (IKC). Take into account shoe heel height, and add 3mm safety factor.

b) Ideally, the pylon connecting the knee and foot should end up vertical. Of course, there may be a variance due to the foot alignment recommendations. In this case, the maximum

! Set the bench alignment taking into account the heel height of associated footwear plus 3mm safety factor!

! It is not recommended to have alignment posterior to the reference line, as it could cause knee instability!

anterior tilt of the pylon should not to exceed 4 degrees, and it may be necessary to utilize offset adapters like the 1222T off set tube clamp. In some cases, it may be necessary to slightly adjust the pyramid angle adjustment by loosening the set screw for the adjustment screw and then tuning the angle adjustment screw clockwise. **⚠ DO NOT OVERLY ADJUST THIS SETTING** – it will reduce the K point or IKC making the knee more unstable. Some adjustment is acceptable, but take into account all aspects of patient ability, length of limb, and foot when going to this adjustment. It is advised to follow up in 1-2 weeks to reassess the alignment.

c) The weight line should pass through the centerline of the knee in the Coronal or M/L plane (Fig. 4). Excessive outset or inset will put undue stress on the knee joint.

d) The weight line for Sagittal or A/P plane should have the plumb line passing through T (Tg) line. Ideally, Tg line should pass through the knee center (red colored pivot) and be perpendicular to the ground. (Fig. 3)

e) **When to use ZONE 1 option** - For the higher weight spectrum patients, Tg line should pass slightly into the “Zone I” area (up to 2mm), which is indicated as up to “T1” (Fig. 3) so that it will reduce stance flexion moment forces of the 5th bar. It is recommended to have controlled stance flexion action of the 5th bar, and not to lock it out or have excessively long duration of it. Excessive stance flexion moment and/or duration can be adjusted through tightening the “Flexion Control Adjusting Screw” in, and/or adjusting alignment (of socket and/or foot) into the T1 zone (up to 2mm). The goal is to have smooth transition from stance flexion to neutral mid-stance motion, and not excessive stance flexion duration.

f) **When to use ZONE 2 option** - For the lighter weight spectrum patients, such as women and children, Tg line should pass slightly into the “Zone II” area (up to 2mm), which is indicated as up to “T2” (Fig.3) so that it will increase stance flexion moment forces of the 5th bar linkage. It is recommended to have controlled stance flexion action of the 5th bar, and not excessively short duration, or abrupt motion of it. Excessive short stance flexion moment and/or duration can be adjusted through loosening the “Flexion Control Adjusting Screw” out, and/or adjusting alignment (of socket and/or foot) into the T2 zone (up to 2mm). The goal is to have smooth transition from stance flexion to neutral mid-stance motion, and not abrupt motion or excessively short stance flexion duration.

CAUTION: Please pay extra caution on Tg line passing towards the maximum (up to 2mm) or past Zone I because it will cause excessive knee head extension force, which will generate excessive leverage pressure on "back linkage" and "fifth axial bar" and could result in knee breakage as shown in Fig. 5.



It is highly recommended that Tg line should ONLY pass through the "Red Dot". (Fig. 3)



Fig. 5 An exaggerated schematic diagram to show affect on 5-bar linkage with excessive anterior alignment

8 Knee Adjustment

8.1 Flexion / Extension Adjustment

Swing Phase Pneumatic setting is pre-set from the factory. Extension or flexion adjustment is only needed if the clinician finds the wearer shows a need for higher or lower walking speeds.

Swing Phase Control Adjustment: It is advisable to adjust **flexion before extension** for optimum walking symmetry. If needed, please follow directions below.


! Ensure full knee extension occurs before performing extension adjustment.

! Use following procedure only if there is a need to adjust extension or flexion:

- 1) Turn extension screw anti-clockwise to lowest resistance - then
- 2) Turn flexion screw clockwise to set to highest resistance

! (Do not over tighten if screw has resistance, or bottoms out – damage may occur and void warranty!)

- Incrementally loosen (anti-clockwise) the flexion screw to adjust heel lift;
- Incrementally tighten (clockwise) the extension screw to smoothly stop extension.

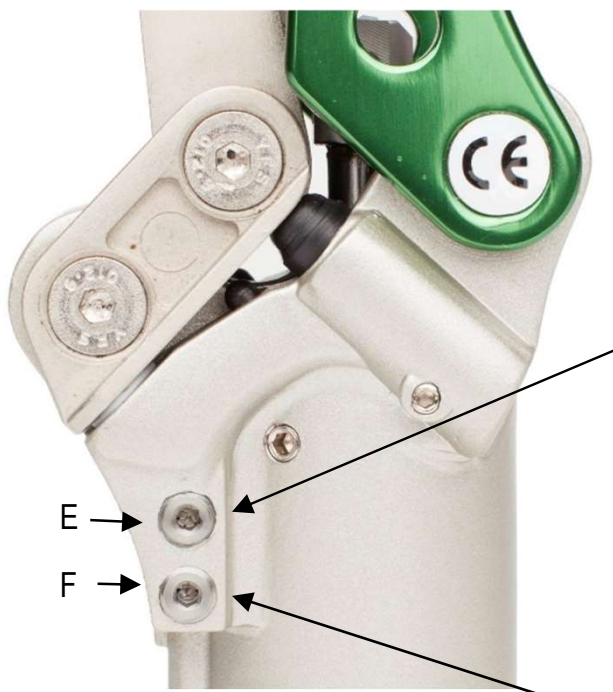

Extension Resistance Adjustment: (Top Screw)

Turn Extension adjustment screw with 2.5mm hex wrench:

Clockwise to increase extension resistance.


Anti-clockwise to reduce extension resistance.

Turn Flexion adjustment screw with 2.5mm hex wrench: (Bottom Screw)

Clockwise to increase flexion resistance.

Anti-clockwise to reduce flexion resistance

8.2 Stance Flexion Adjustment

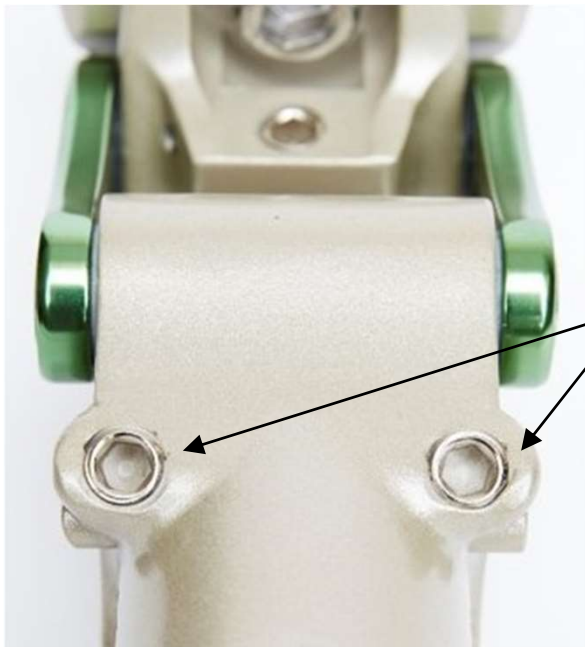

NOTE: Stance Flexion set screws have been eliminated. The threaded holes remain, but will not contain a screw.




5mm Stance flexion adjustment screws are located on the anterior body of the knee. Turning both screws clockwise decreases the 5th axis motion, reducing the stance flexion angle. Anti-clockwise will increase 5th axis motion, increasing the stance flexion angle. Adjustments need to be symmetrical!




Important! adjustments to Stance Flexion Bumper screws need to be symmetrical! Apply thread locker to screws one at a time to prevent screw backing out.

8.3 Pyramid Head Position Adjusting

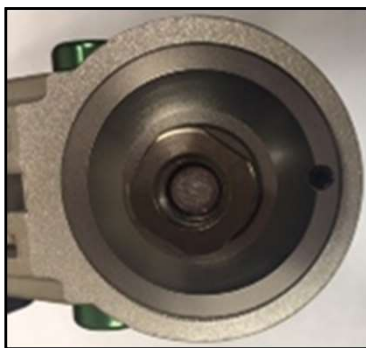


Loosen screw with 2.5mm hex wrench.

Loosen Pyramid bolt with 8mm hex wrench. Rotate to desired orientation And retighten bolt.

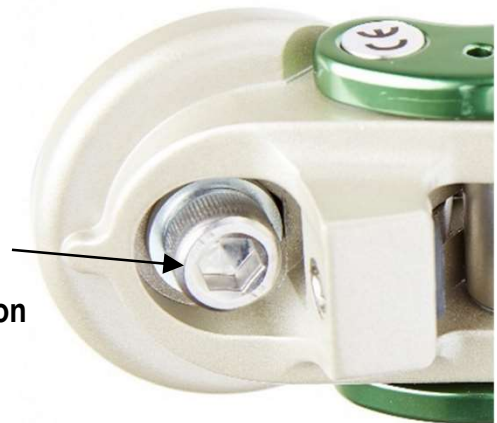

Note: Mark/indicate pyramid orientation. Remove pyramid bolt and apply thread locker, and torque bolt 18Nm. Tighten set screw to help prevent rotation.

8.5 Knee Head Tilting Adjusting

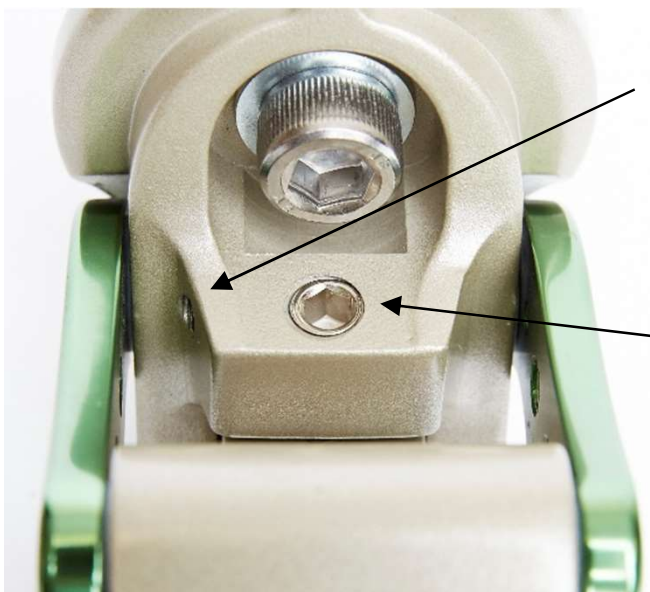

Loosen the set screw prior to adjusting. After adjusting, tighten set screw.

Use 5mm driver and turn the screw:
Clockwise (In) to reduce geometric stability of the knee (high front/low rear) and will make knee flexion easier.
Anti-clockwise (Out) to increase geometric stability of the knee (low front/high rear) and will make knee flexion more difficult.



Level or very slight up-sloping tilt (Factory Setting) is most common. Excessive adjustment from factory setting will cause decrease in bumper life span due to excessive compression forces from adjustment screw and decreased knee stability!

Pyramid Head Tilt Adjustment:




Adjusting the Head Tilt will move the IKC. If the IKC is set too anterior (less stable) and the socket weight line too far posterior, the knee will be very unstable and could lead to fall or injury. Head Tilt is to “Fine Tune” the knee and will not substitute for inadequate socket placement!

Anti-clockwise rotation will tilt the pyramid anteriorly which increases stability by adding extension into the socket and moving the weight line anteriorly versus the ICR. (ICR is more posterior of socket making flexion initiation harder)


Clockwise rotation will tilt the pyramid posteriorly which decreases stability by adding flexion into the socket and moving the weight line posteriorly versus the ICR. (ICR is more anterior of socket making flexion initiation easier)

IKC

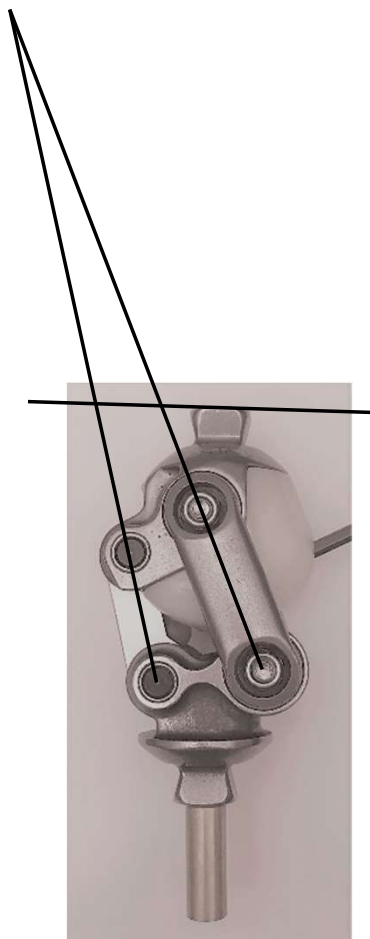

Very Stable
Harder to initiate swing

IKC

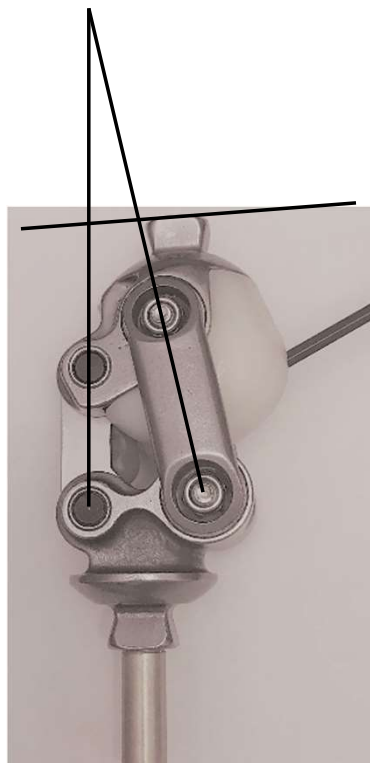

Stable
Neutral Function

IKC



Less Stable
Easier to initiate swing

9. Knee Maintenance

9.1 Changing Various bumpers:

For the knee head leveling bumper, use 2mm hex wrench to loosen the bumper set screw by turning anti-clockwise. Fig. right



Note: Indicate position of the screw prior to removal. Push out the bumper by turning the screw with 5mm hex wrench clockwise. Photo left

Place new bumper back into the slot same as previous one.

For the round shape bumper, use a small standard screwdriver to pry out the old bumper. Insert with a new one.



9.2.2 Changing Stance Flexion Bumper

Use 4mm hex wrench to loosen the screws of fourth and fifth axes to take out the side bars of stance flexion unit. Please refer to the picture below




Pull out opposite side bar along with fourth and fifth axes. Stance flexion bumpers should be accessible for servicing.



10 Technical Specification

- | | |
|--|--|
| · Operating & Storage Temperature Range: | -10°C to 50°C (14°F to 122°F) |
| · Weight: | Pyramid / Lotus 877g / 892g |
| · Recommended Activity: | K2 |
| · Maximum User Weight: | 125kg (275lbs) |
| · Maximum flexion angle: | 135 degrees |
| · Proximal Alignment attachment: | Rotatable Male Pyramid, or Lotus adapter |
| · Distal Alignment attachment: | Tube Clamp |
| · Tube clamp torque setting: | 12Nm |
| · Pyramid Center Bolt: | 18Nm |
| · Build: | Pyramid / Lotus 177.5 / 183.8 |
- **Materials: Aluminum Alloy, Stainless Steel, Steel, Rubber**
-

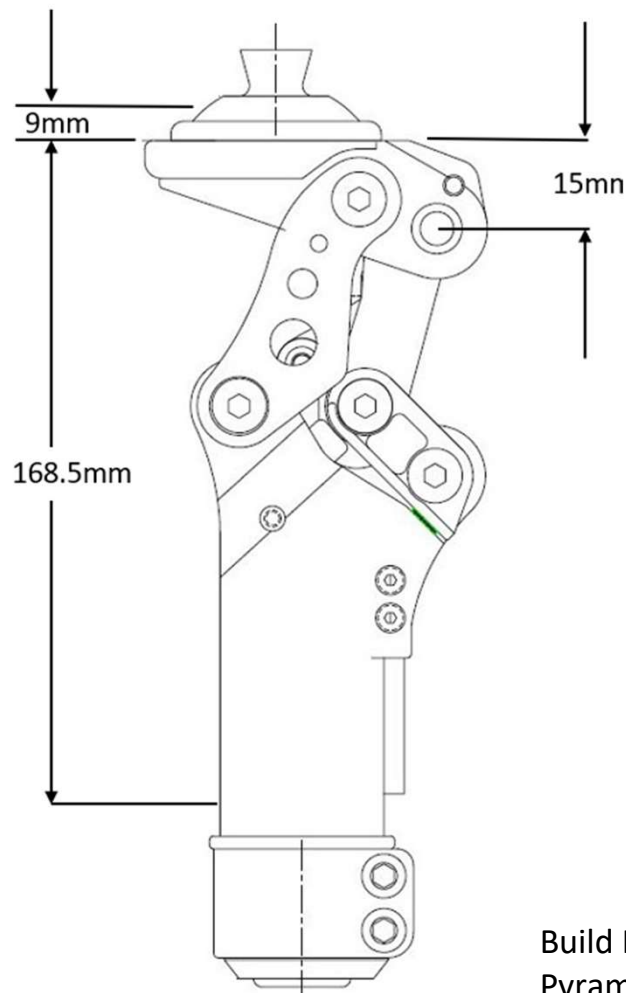

Build Height:
Pyramid = 177.5mm
Lotus = 183.8mm

10 Warranty

Warranted for 2 years from the date of invoice by ST&G.

The user should be aware that changes or modifications not approved will void the warranty.

11 Liability

The manufacturer recommends using the device only under the specified conditions and for the intended purposes. The device must be maintained according to the instructions for use supplied with the device. The manufacturer is not liable for damage caused by the component combinations that were not authorized by the manufacturer.



CE Conformity

This product meets the requirements Council EU 2017/745 (EU) for medical products. This product has been classified as a class I product according to the classification criteria outlined in appendix IX of the guidelines. Please keep this manual in safe place for future use.

Useful Life

Service life of the product is covered by the warranty period.

Disposal

The device and its packaging must be disposed of in accordance with respective national/local environmental regulations.

1323 IFU REV F (4-30-25)