



4-Hole Rotating Base

Instructions For Use

PA13




1. Application & Description

- This adapter is exclusively for lower extremity prosthetic fittings
- It connects any 13 series knee with mounting platform to adapt to a pyramid receiver
- Compatible with:
 - 1321, 1322, 1323, 1324, 1326, L1321, L1324

2. Technical Specifications

P/N	Material	Build Height	Weight Limit
PA13	Aluminum / Stainless Steel	10mm	125 kg / 275 lbs

- **8 mm Hex Bolt Torque:**
 - 30 Nm
- **2.5 mm Hex Rotation Set Screw Torque:**
 - 7 Nm
- **Thread Locker:** Apply Blue Loctite

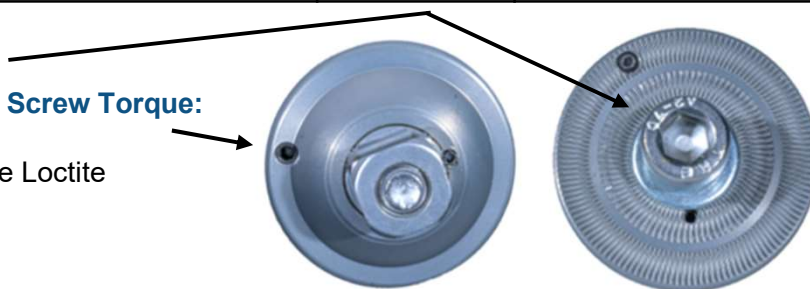

3. Installation Instructions

Application

- **Trial Alignment:** Place the PA13 onto the serrated knee mounting platform. Adjust alignment and secure the mounting bolt with an 8 mm hex wrench torqued to 30 Nm for dynamic alignment.
- **Mark Position:** Once alignment is finalized, trace the outline of the PA13 base and mark the rotation set screw location on the mounting platform.
- **Secure Mounting Bolt:** Remove the mounting bolt, apply blue Loctite, re-install, and torque to 30 Nm.
- **Lock Rotation:** Apply blue Loctite to the rotation set screw and torque to 7 Nm using a 2.5 mm hex wrench.

4. Safety & Environment

- **Temperature Range:** Between -20°C to 80°C (-4°F to 176°F).
- Rinse with clean water if exposed to salty, chlorinated, or otherwise potentially corrosive liquids.
- **Storage:** Store at room temperature; avoid wet, very hot, or extremely cold environments
- **Warranty:** 2 years, provided it is used as intended without modification, and only for the original intended user.
- Do not reuse on another prosthesis.
- **Conformity:** This product meets the requirements Council EU 2017/745 (MDR) for medical products. This product has been classified as a class I product according to the classification criteria outlined in appendix IX of the guidelines. Please keep this manual in safe place for future reference.



PA13 IFU REV A (9-23-2026)
Specifications subject to change.



MDSS GmbH