

# GO MAGNETIC MAGLOCK 2.0 & 2.0 PLUS

INTRODUCING THE NEWLY DESIGNED  
LOW PROFILE, NO BATTERY – NO WIRES

The **MAGLOCK & MAGLOCK PLUS** have been designed to ease the frustration of donning a pin lock liner for patients K1 to K3. The **MAGLOCK 2.0** achieves up to 60lbs of suspension strength while the **MAGLOCK 2.0 Plus** achieves up to 90lbs. The large distal target allows the patient to be 10 degrees off center and still achieve suspension.

Neither version requires any maintenance, no charging and offers a fail safe way to turn on and off the magnetic system. The ultimate product for patients who have dexterity issues, issues with pin alignment, requires a caregiver or looking for a suspension system that offers an easier option for suspension.



\* Patented



MAGLOCK 2.0 PLUS

MAGLOCK 2.0




Push and slide mechanism

P/N	Product Description	Weight Limit	Product Weight	Build Height	Suspension Strength
3180M	Mechanical MAGLOCK	136 kg / 300 lbs	278g / 0.6 lbs	2.2 cm / 0.87 in	60 lbs
3181MP	Mechanical MAGLOCK 2.0 PLUS	136 kg / 300 lbs	312 gm / 0.69 lbs	2.72cm / 1.07 in	90 lbs
3180-1M	DUMMY for Maglock 2.0				
3181-1M	DUMMY for MAGLOCK 2.0 PLUS				
DSTC	Distal Cap Converter				

## DISTAL CAP Converter for all pin liners

Cap converter with 10 mm thread for all pin liners. The distal cap converter is required for use with other liners in order to be compatible with the MAGLOCK. DSTC comes with two silicone rings

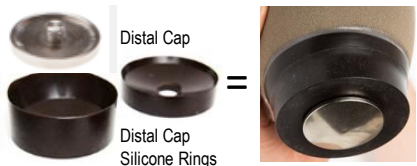

## Fabrication Dummy for MAGLOCK & MAGLOCK PLUS


Comes with different sized silicon rings for casting over a MAGLiner or when using the Distal Cap converter for a locking liner. These rings provide the perfect spacing and alignment of the MAGLOCK. The Fabrication Dummy comes with two silicone rings.



# TRANSTIBAL MAGLiners

## PREMIER & HYDRO

### MAGLiner - Premier

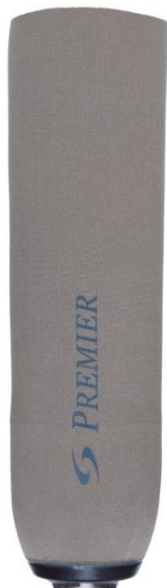

This premium liner is designed for the health of the residual limb and ease of donning. Manufactured with mineral oil embedded Platinum cured medical grade silicone and sheathed in a seamless durable fabric. The matrix which is built into the liner prevents distal elongation.

6 Month warranty

MAGLiner – PREMIER			
P/N	Description	P/N	Description
3040-18	MAGLiner Premier, 18	3040-28	MAGLiner Premier, 28
3040-20	MAGLiner Premier, 20	3040-30	MAGLiner Premier, 30
3040-22	MAGLiner Premier, 22	3040-32	MAGLiner Premier, 32
3040-23.5	MAGLiner Premier, 23.5	3040-34	MAGLiner Premier, 34
3040-25	MAGLiner Premier, 25	3040-36	MAGLiner Premier, 36
3040-26.5	MAGLiner Premier, 26.5		

\* Using the MAGLINER with the MAGLOCK offers the lowest profile possible.

### MAGLiner - Hydro



The Hydro MAGLiner's silky smooth Platinum cured medical grade silicone liner has the unique ability to repel water. This hydrophobic property makes the Hydro the easiest and fastest liner to clean and sanitize. The matrix which is built into the liner prevents distal elongation. The Parylene® surface reduces surface friction between the liner/skin and liner/socket for better comfort and donning. The Hydro MAGLiner reduces skin shear for fragile skin and will help with blister reduction. The Hydro Skin offers proprietary donning and doffing strips to aid the patient in easier application/removal. It is suitable for moderate to high activity patients who require cushioning over sensitive areas in the socket. The Hydro Skin liner's opaque pigment gives a natural "skin tone" to all ethnic skin types.

6 Month warranty

MAGLiner – HYDRO			
P/N	Description	P/N	Description
3030-18	MAGLiner Hydro, 18	3030-28	MAGLiner Hydro, 28
3030-20	MAGLiner Hydro, 20	3030-30	MAGLiner Hydro, 30
3030-22	MAGLiner Hydro, 22	3030-32	MAGLiner Hydro, 32
3030-23.5	MAGLiner Hydro, 23.5	3030-34	MAGLiner Hydro, 34
3030-25	MAGLiner Hydro, 25	3030-36	MAGLiner Hydro, 36
3030-26.5	MAGLiner Hydro, 26.5		

\* Using the MAGLINER with the MAGLOCK offers the lowest profile possible.

**MAGLiners, PREMIER & HYDRO are designed exclusively for the MAGLOCK Suspension System.** The MAGLOCK eases the frustration of donning a pin lock liner for patients K1 to K3. The MAGLOCK is developed using a rare earth magnet which achieves up to 90lbs of adhesion. The large distal target allows the patient to be 10 degrees off center and still get suspension.



**ST&G CORPORATION**  
2691 Saturn Street  
Brea CA 92821, U.S.A.

TO ORDER CALL: 714-524-0663 / FAX: 714-364-8113

# TRANSFEMORAL MAGLiners

## PREMIER & HYDRO

### MAGLiner - Premier



This TF premium liner is designed for the health of the residual limb and ease of donning. Manufactured with mineral oil embedded Platinum cured medical grade silicone and sheathed in a seamless durable fabric. The matrix which is built into the liner prevents distal elongation.

6 Month warranty

MAGLiner – PREMIER			
P/N	Description	P/N	Description
3060 TF - 30	MAGLiner Premier, 30	3060 TF - 38	MAGLiner Premier, 38
3060 TF - 32	MAGLiner Premier, 32	3060 TF - 40	MAGLiner Premier, 40
3060 TF - 34	MAGLiner Premier, 34	3060 TF - 45	MAGLiner Premier, 45
3060 TF - 36	MAGLiner Premier, 36		

\* Using the MAGLINER with the MAGLOCK offers the lowest profile possible.

### MAGLiner - Hydro



The TF Hydro MAGLiner's silky smooth Platinum cured medical grade silicone liner has the unique ability to repel water. This hydrophobic property makes the Hydro the easiest and fastest liner to clean and sanitize. The matrix which is built into the liner prevents distal elongation. The Parylene® surface reduces surface friction between the liner/skin and liner/socket for better comfort and donning. The Hydro MAGLiner reduces skin shear for fragile skin and will help with blister reduction. The Hydro Skin offers proprietary donning and doffing strips to aid the patient in easier application/removal. It is suitable for moderate to high activity patients who require cushioning over sensitive areas in the socket. The Hydro Skin liner's opaque pigment gives a natural "skin tone" to all ethnic skin types.

6 Month warranty

MAGLiner – HYDRO			
P/N	Description	P/N	Description
3050 TF - 30	MAGLiner Hydro, 30	3050 TF - 38	MAGLiner Hydro, 38
3050 TF - 32	MAGLiner Hydro, 32	3050 TF - 40	MAGLiner Hydro, 40
3050 TF - 34	MAGLiner Hydro, 34	3050 TF - 45	MAGLiner Hydro, 45
3050 TF - 36	MAGLiner Hydro, 36		

\* Using the MAGLINER with the MAGLOCK offers the lowest profile possible.

**MAGLiners, PREMIER & HYDRO are designed exclusively for the MAGLOCK Suspension System.** The **MAGLOCK** eases the frustration of donning a pin lock liner for patients K1 to K3. The **MAGLOCK** is developed using a rare earth magnet which achieves up to 90lbs of adhesion. The large distal target allows the patient to be 10 degrees off center and still get suspension.



**ST&G CORPORATION**  
2691 Saturn Street  
Brea CA 92821, U.S.A

TO ORDER CALL: 714-524-0663 / FAX: 714-364-8113